


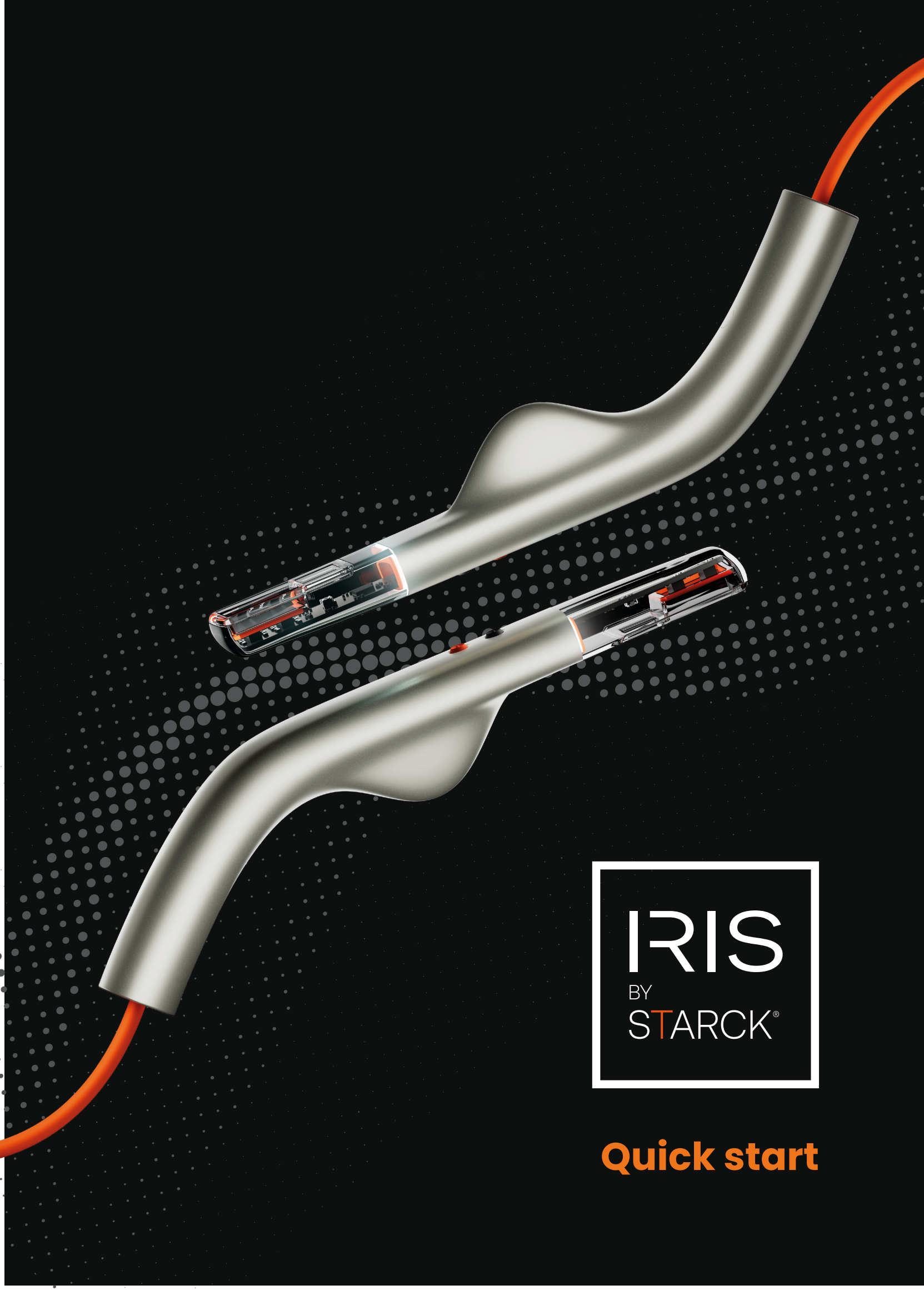


IRIS
BY
STARCK®

Quick start

CONGRATULATIONS !

You have just purchased IRIS BY STARCK®, a comprehensive digital impression solution designed for oral health professionals. This medical device combines an intraoral scanner with powerful proprietary software to create a 3D model of the entire oral cavity.

Enjoy a secure, open, scalable, and interoperable platform thanks to standard formats (including STL and PLY) that allow you to centralize the entire dental treatment workflow: from initial examination to prosthetic fabrication, including patient file management, collaborative work, inter-professional exchanges and order tracking.

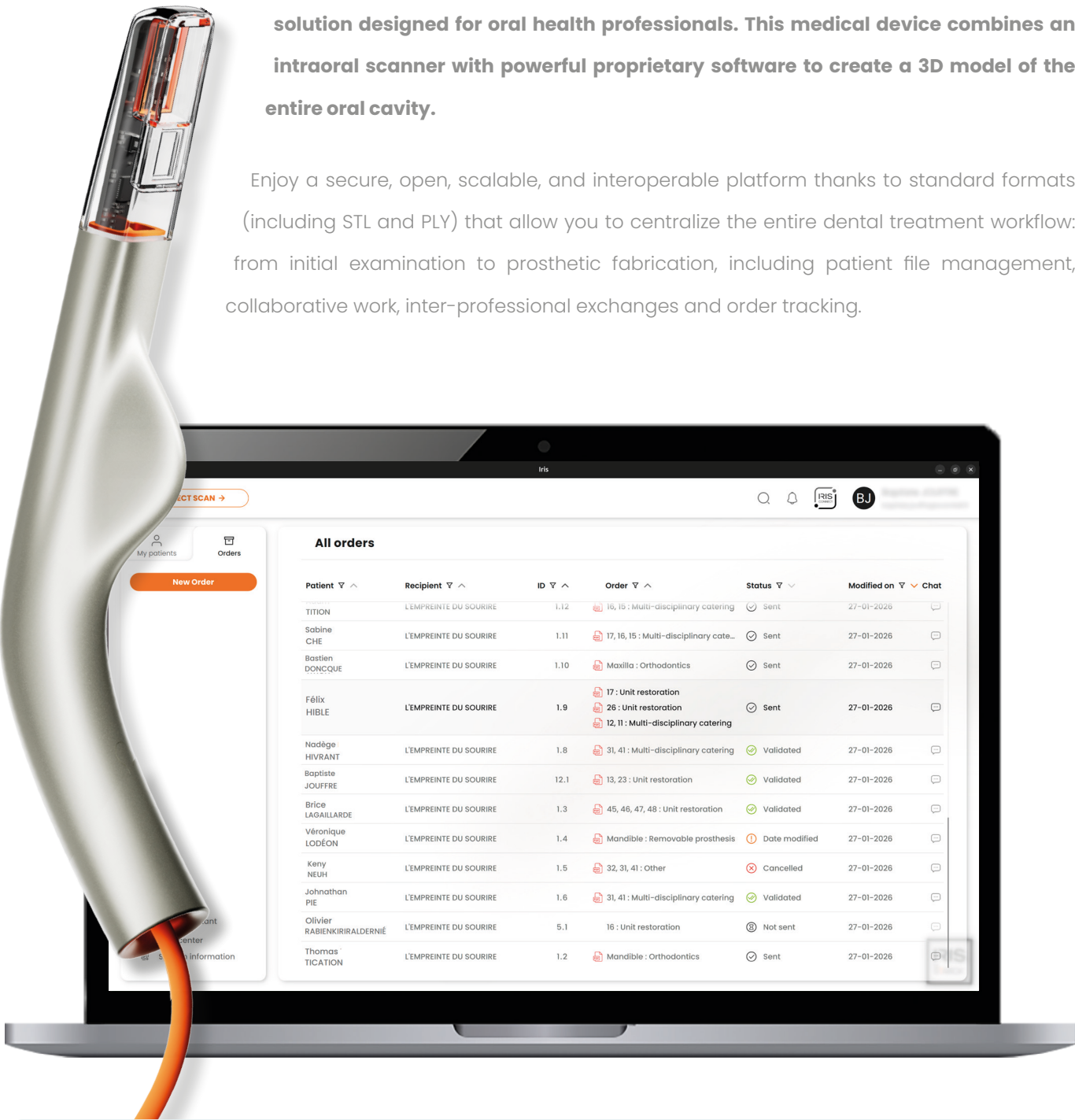

So why "BY STARCK"?

Developed and manufactured by Scan4all, the handpiece was designed by Philippe Starck. It stands out for its organic, biomimetic lines and a shape adapted to the user's morphology.

CONTENTS

- DEVICE OVERVIEW 4
- ACTIVATE YOUR LICENCE 4
- YOUR HOME PAGE 5
- A CUSTOMIZED JOURNEY 6
- REAL TIME SCAN 7
- TOOLS RELATED TO THE SCAN 8
- MY PATIENTS 10
- MY ORDERS 11
- IRIS CONNECT 12
- MY ACCOUNT 14
- MY SETTINGS 15

DEVICE OVERVIEW



1 x ready-to-use laptop



1 x Handpiece IRIS BY STARCK®



1 x SECURE.ID card with your login credentials



1 x Support base



1 x USB flash drive



4 x Pairs of black glove



1 x Microfiber cloth



1 x Disinfection base

ACTIVATE YOUR LICENCE

Have you received a computer installed by SCAN4ALL?

Launch IRIS BY STARCK®, click on "activate my account" and enter your license number when registering using the unique number printed on the SECURE.ID card.





Plug the computer into the mains.



Connect the handpiece to a USB-C port.



Turn on the computer.
Use the password written on your card included in the package.



Connect to the internet.




Launch IRIS BY STARCK® and activate your license when you register.



Can I use IRIS BY STARCK® without an internet connection? Yes, once activated, you can scan and use the IRIS software without any restriction. Your orders will be sent automatically once an internet connection is established.

YOUR HOME PAGE

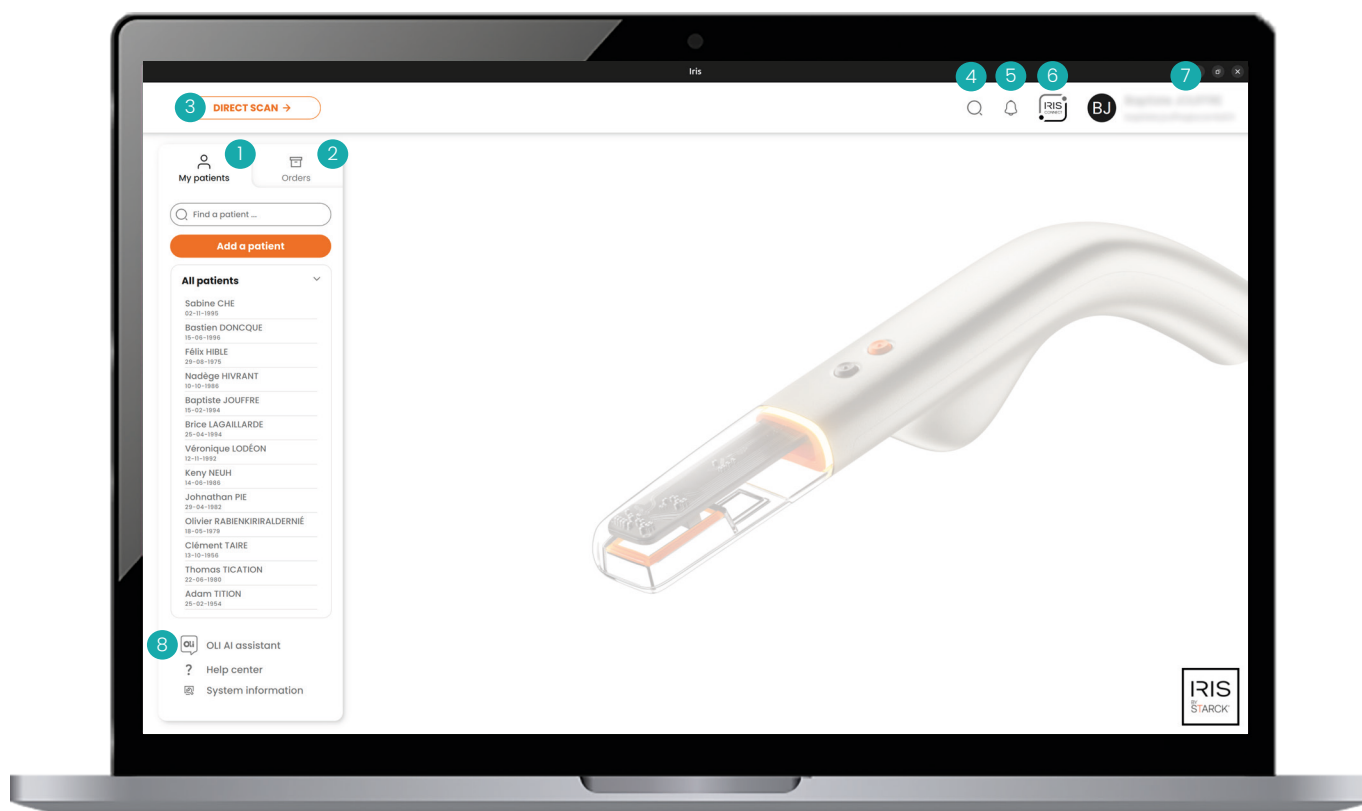

Enjoy a completely free, autonomous, and user-friendly ecosystem that naturally supports you in your daily life.

- 1 **Manage the list of all my patients** : search and add.
- 2 **Manage my orders** : creation, tracking and history.
- 3 **Start a scan** without prior information.
- 4 **Search bar** (patients, scans, orders, contacts, collaborators, or assistants).
- 5 **Notification center** for actions related to your contacts and orders.
- 6 **IRIS CONNECT** : collaborate, discuss, share, manage my contacts and my team.
- 7 **My account and settings**.
- 8 **Chat with OLI**, an integrated artificial intelligence that provides level 1 technical support.

A CUSTOMIZED JOURNEY

Whether you start with a scan or an order, the user interface adapts automatically to your logic. Switch between workflows at any time without losing information.



The workflow corresponds to the different steps required to send an order : from creation of the prosthetic file to the sending of the selected scans.

For those who prefer to scan first, then organize

DIRECT SCAN →



Scan your patient immediately using the **DIRECT SCAN** feature. Then, save the scan to a patient record and create an order later.

DIRECT SCAN is a flexible workflow. This mode prioritizes responsiveness while retaining all standard features.

For those who plan before acting



Edit the prosthetic file > Maxillary scan > Mandibular scan > Occlusion scan > Send order.

The **GUIDED workflow** automatically guides you through the process based on the treatment selected in the prosthetic record.

Each type of procedure **triggers the scanning steps required for the order** : all you have to do is follow the workflow.

You retain the freedom to customize this workflow.

The **OPEN workflow**, where each step is independent (with a click).

You are responsible for performing the scans necessary for a complete order by selecting the steps yourself.



From the "Orders" tab or from a patient file, click on "New Order." Fill out your prosthetic form and choose your preferred workflow mode (guided or open).

REAL TIME SCAN



Connect the computer to a power outlet.



Avoid direct light on the area to be scanned.



Use black gloves.



Dry the mouth and teeth.



Start scanning.



To begin your scan, you can start with the occlusal surface, then move on to the lingual or palatal surface, and finish with the buccal surface.



For incisors, rock back and forth on the occlusal axis and finish by sliding slightly on the front surface.

Scan at the correct distance: automatic calibration

Positioning the tip parallel to and in contact with the teeth ensures the correct distance for an accurate scan.

Scan at a steady pace: maintain a fluid motion

Adjust your speed to match the refresh rate of the 3D view and prioritize a steady motion, *like a paintbrush gliding across a canvas*. The IRIS BY STARCK® scanner *glides over the teeth smoothly*, without jerks or rushing.

Light indicator

On or off: white, the camera is on.

Tracking status: same as the scan view.

Green : tracking in progress.

Yellow : slow down; the scan is stabilizing.

Red : loss of tracking.

Start the scan.

Pause.

Resume.

Confirm and start the final reconstruction.



To ensure a smooth scanning experience: use the timer icon (in the top right corner of your real-time scan view).

Gray: You can continue.

Yellow: Scanning experience is slowed down.


Red: Prioritize starting the final reconstruction. You can resume scanning using the “Resume Scan” tool.

TOOLS RELATED TO THE SCAN

IRIS BY STARCK® offers a suite of tools to control, adjust, analyze, and enhance your scans.

Refine your visualization

Take advantage of the **visualization features** 1 to view PLY or STL files, take screenshots, change your scan viewpoint, adjust transparency, and choose which measurements to display or hide.

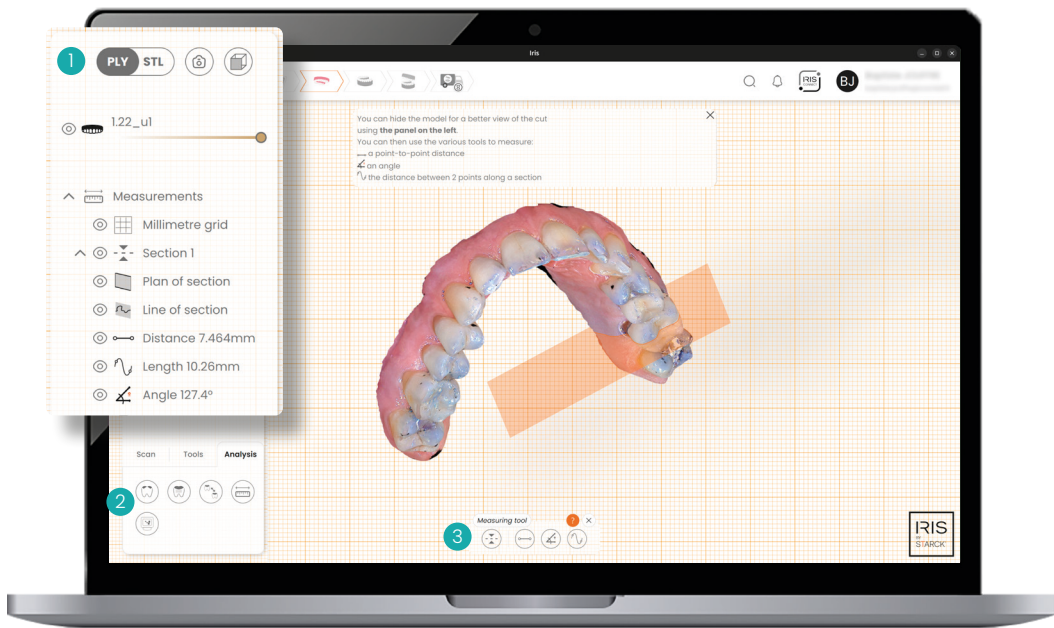


Access **advanced analysis, editing, and customization functions** 2 through three complementary tabs : Scan, Tools, and Analysis.

 Draw a section 3 on a 2D plan to make it easier to locate and display the millimeter grid for even greater precision.

Control every detail of your scan

Take advantage of **scan management features**, **tools** to support your clinic, and capabilities for comparative analysis and measurement.

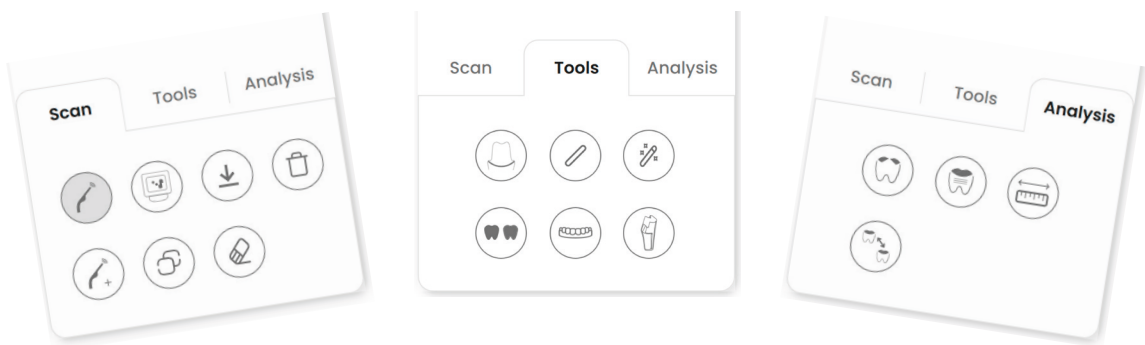

All available tools

SCAN

 Resume scan

 Add a scan

 Copy scan

 Erasing the scan

 Download scan

 Delete scan

 Photogrammetry

TOOLS

 Display finishing lines

 Create a manual finishing line

 Create an auto finishing line

 Change the type of base

 Label

 Add a scanbody

ANALYSIS

 Curvature

 Roughness

 Comparison

 Measures

Even greater precision in occlusion

1

Measurements

Millimetre grid

Section 1

Plan of section

Line of section

2

Analysis

5.1_ul

5.1_ll

3

Color scale

1.00

0.50

0.15

0.05

Rugosité

 Add a scan to occlusion

 Optimisation of intercuspidity

 Add an occlusion

 Adjust manually

Find all the **visualization features** 1 and previous analysis tools, as well as **additional features specific to occlusion** 2.


Optimisation of intercuspidity

Automatically perfect the intercuspidity thanks to alignment based on previously identified contact areas.


3 Customize the value scale

Color codes help you quickly locate critical areas. You can define your own colorimetry.

MY PATIENTS

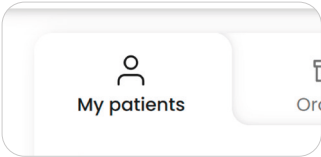

On the home page, the “My Patients” tab allows you to quickly access a list of all your patients, add new ones, view and manage their information, or view their treatment history (see below).

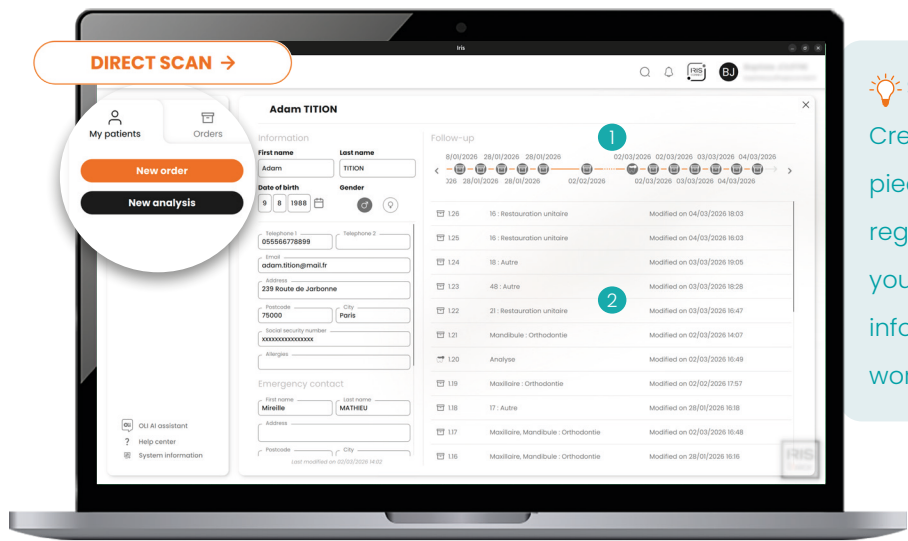

Save time

Creating a patient file requires only 4 pieces of information, allowing you to register a new patient in seconds, and your assistants can fill in additional information without interrupting your workflow.

By selecting a patient, you can view an optimized overview of their treatment history on their timeline 1 and in a chronological list 2. Each action is recorded in the patient file.

3 types of procedures

DIRECT SCAN

ORDER

ANALYSIS

Stay in control

Rights	None	Reading	Writing
Patient file	☐	☒	☐
Discussion	☒	☐	☐

CancelConfirm

Set up custom access rights to allow your assistants to fill in the information (see “IRIS CONNECT”, page <CM> 12).

Quick search


Search ...

From the search bar in the tab or via the general search bar: quickly locate your patient to access their medical history.

Permanent deletion : When a patient is deleted, all of their information and associated scans are permanently deleted. No records are retained in accordance with applicable legislation.

MY ORDERS

On the home page, the "Orders" tab lets you **track all your orders**. You can place an order at any time directly from this space.

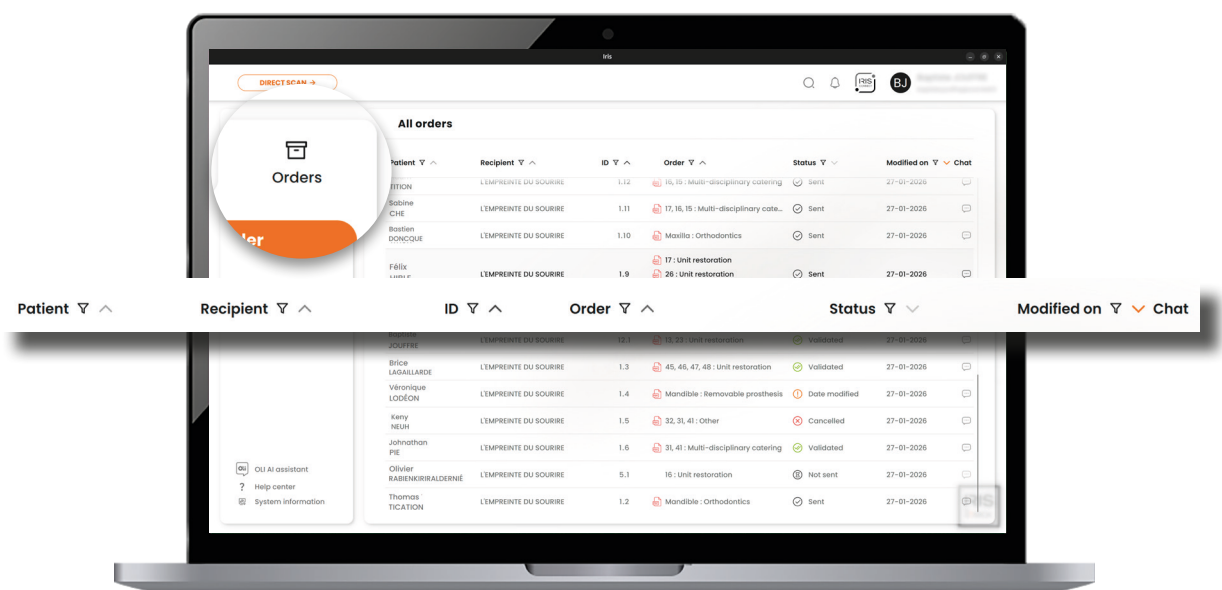

💬 A dedicated chat room to track the discussion for each order.

Place an order

Nothing could be simpler.

From the "Orders" tab or directly via the patient file : click on **"New order,"** fill out the prosthetic form, then select the workflow of your choice. Follow all the necessary steps until the order is sent.

Selecting your partners

Save and send your orders to your partners directly from IRIS BY STARCK® when creating an order :

Synchronization LaGalaxy

Synchronize your data via LaGalaxy account to ensure that your data is continuously updated (enable this in **"My account"**).


Integrated Partners

Choose one of the partners integrated with IRIS BY STARCK® (e.g., SMILERS, CIRCLE) to send the order to your account with that partner.

IRIS CONNECT

Invite your own partners via IRIS CONNECT. They will be able to receive your orders and interact with you via messages (see **"IRIS CONNECT"**, page 12).

IRIS CONNECT

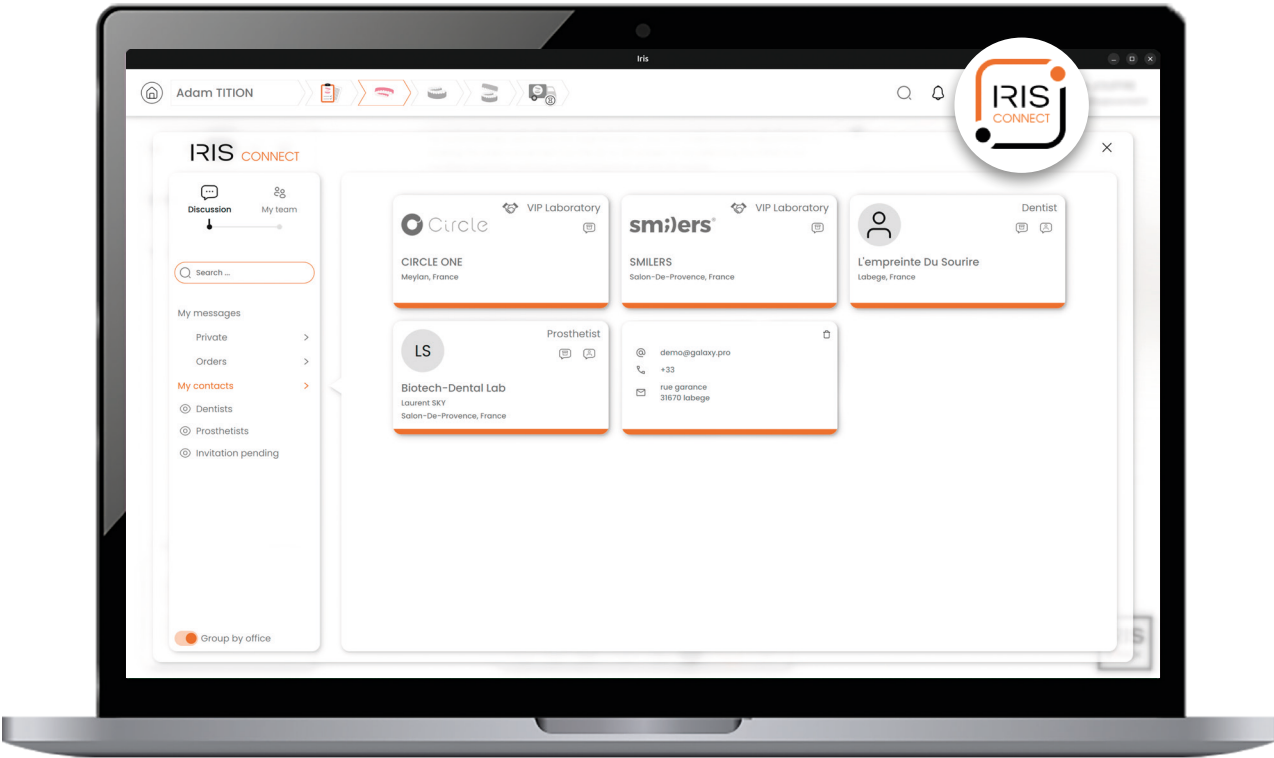

Maximum security

- › Secure sharing environment.
- › ISO 27001 and HDS certification.

Custom team creation

- › Create your team as an administrator.
- › Manage roles tailored to your organization; access rights to patient files and discussions:
 - for your collaborators.
 - for your assistants.

Instant messaging

- › Automatic message management.
- › Smart notifications.



Dedicated discussions forum:

- for your orders.
- for private discussions.

Easy order management


- › Automatic forward of offline orders.
- › Order status.
- › Status update notifications.
- › History of exchanges.



Unlimited collaboration

- › Exclusive IRIS BY STARCK® network directory.
- › Integration of your external contacts.
- › Universal file format (STL).
- › Secure external web interface.

A web interface for your external partners

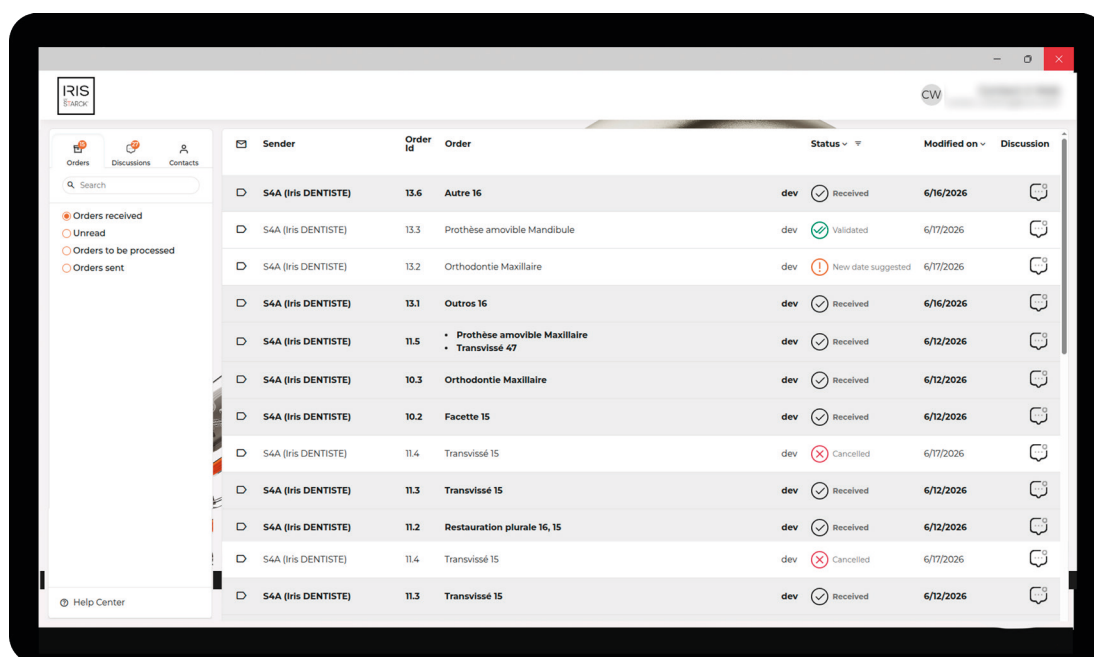

IRIS BY STARCK® WEB is **a universal and secure exchange interface**, free for you and your external partners.

 In addition to having a latest-generation scanner, you can invite any colleague, clinic, or laboratory to discuss the order of your choice.

Invite an external partner

Can't find your contact?

Send an invitation by email

From IRIS CONNECT, use the search field to quickly find a colleague. If they are not listed in the IRIS BY STARCK® directory, **enter their email address to invite them**. They will receive a registration link with a unique security code to create their web account at www.irisbystarck.com

Integrated Partners


IRIS BY STARCK® offers integrated partners: when creating an order, select the partner of your choice as the recipient to view the order directly in your partner account.

Stay informed

Notifications, which you can view by the icon at the top of the home page , keep you updated on activities related to your contacts and orders : reply or send a contact invitation, or be notified of a change in order status (confirmed, canceled, modified).

MY ACCOUNT


The “My Account” tab allows you to **view and modify your personal information** at any time. Convenient and quick, you remain in complete control of your data from your IRIS BY STARCK®.

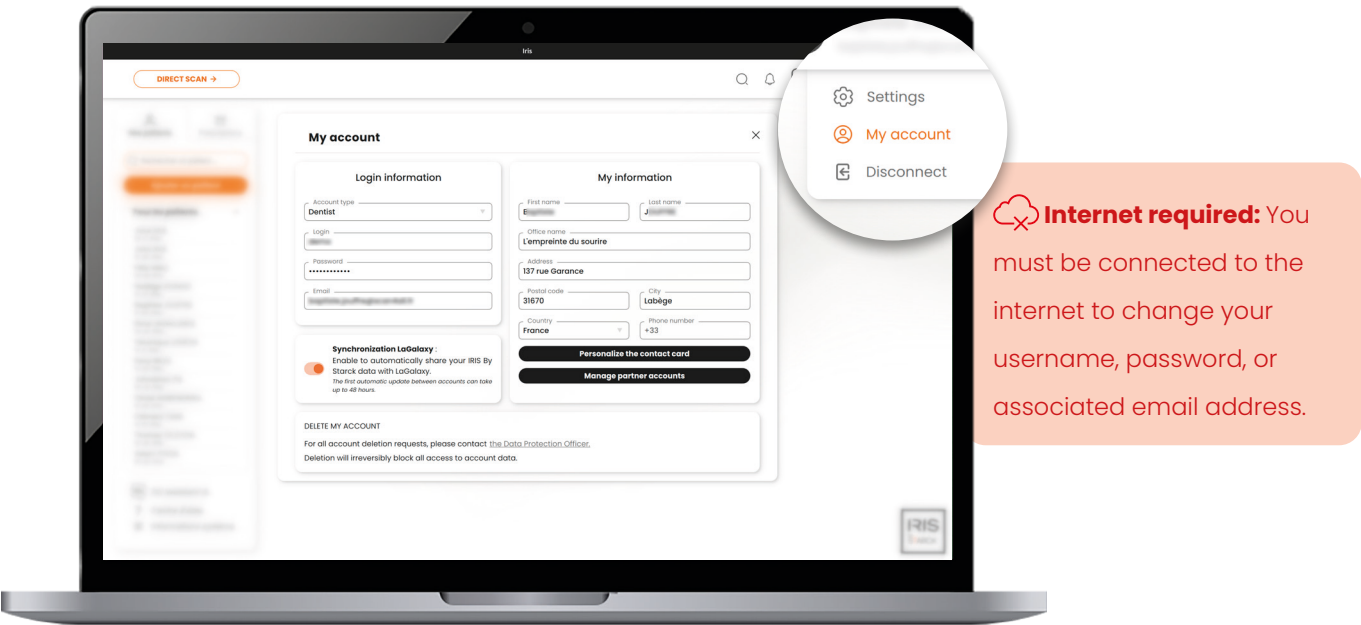

Internet required: You must be connected to the internet to change your username, password, or associated email address.

Your exclusive professional network

Your **visibility on the exclusive IRIS BY STARCK® network belongs to you**. Enjoy the collaborative benefits of IRIS BY STARCK® by joining the community directory, while easily maintaining control over your information.

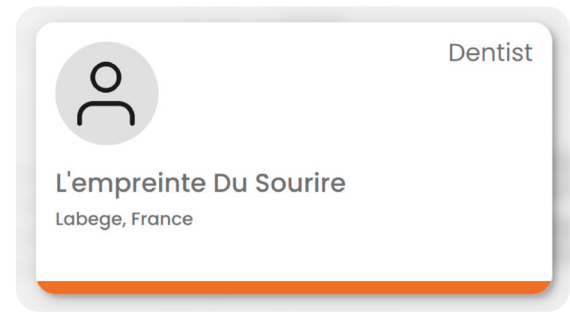


Customize the information displayed on your contact card (automatically generated) : you are free to share or hide the name of your practice, as well as your first and last names.

Connexion with LaGalaxy

Your entire dental world, synchronized on a secure digital ecosystem that brings together all of the Biotech Dental Group’s products and solutions. Once synchronization is activated, your patients, scans, and IRIS BY STARCK® orders will be available in your LaGalaxy space.

MY SETTINGS

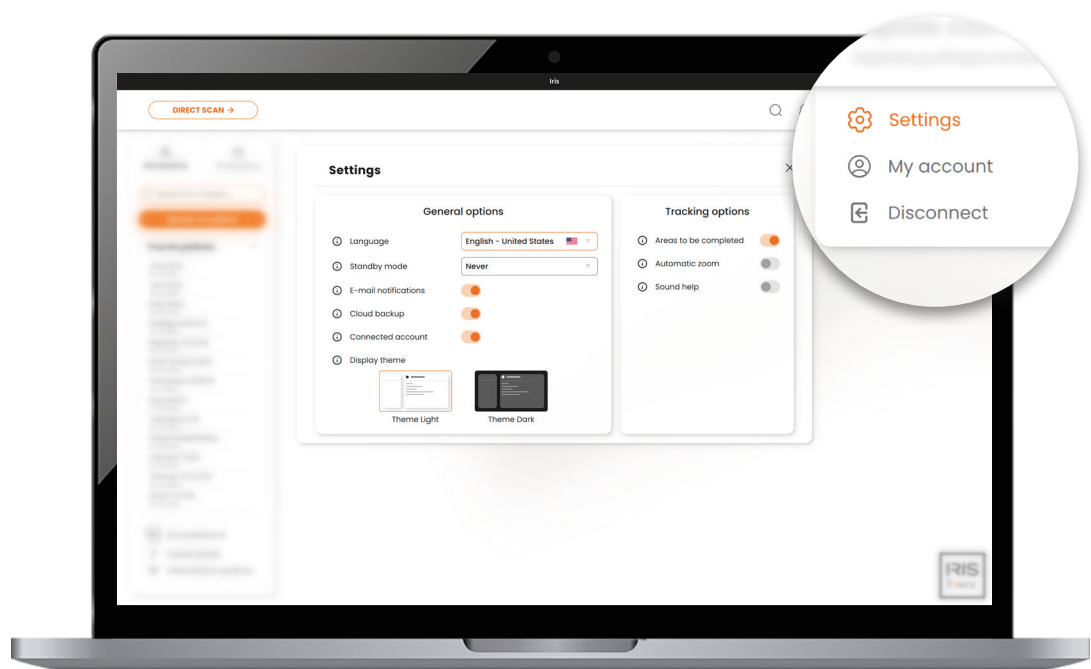

Choice of language, favorite theme, and notifications

Choose the language in which you want to use IRIS BY STARCK®. Adapt your display theme to your environment (light or dark) to ensure maximum visual comfort and choose whether or not to enable email alerts to track order status changes at your own pace.

Automatic standby

Set an inactivity timeout before the application automatically logs out and ask you to re-enter your password to access it.



If you select "Never," remember to log out manually if you are away.

CLOUD backup : free and unlimited

Activate and back up your patient data for free on the IRIS BY STARCK® cloud dedicated to your account, protected by HDS and ISO 27001 security. You must be connected to the internet to use this feature. If you disable it, data will only be backed up on your computer.

Connected account

Use your license on all computers where IRIS BY STARCK® is installed. Enable this feature to retrieve the saved data from all your sessions.



Access all the
documentation

Manufacturer

Scan 4 all

Excellence, for everyone.

137 rue Garance 31670 Labège - FRANCE

N° SIRET 507 484 863 00059

 www.linkedin.com/company/scan4all/

IRIS BY STARCK®

Class 1 medical device (MD) intended for digital dentistry. CE.

For healthcare professionals only.

© Scan4all - Version 1.0 - 2026