


Scan 4 all
introduces

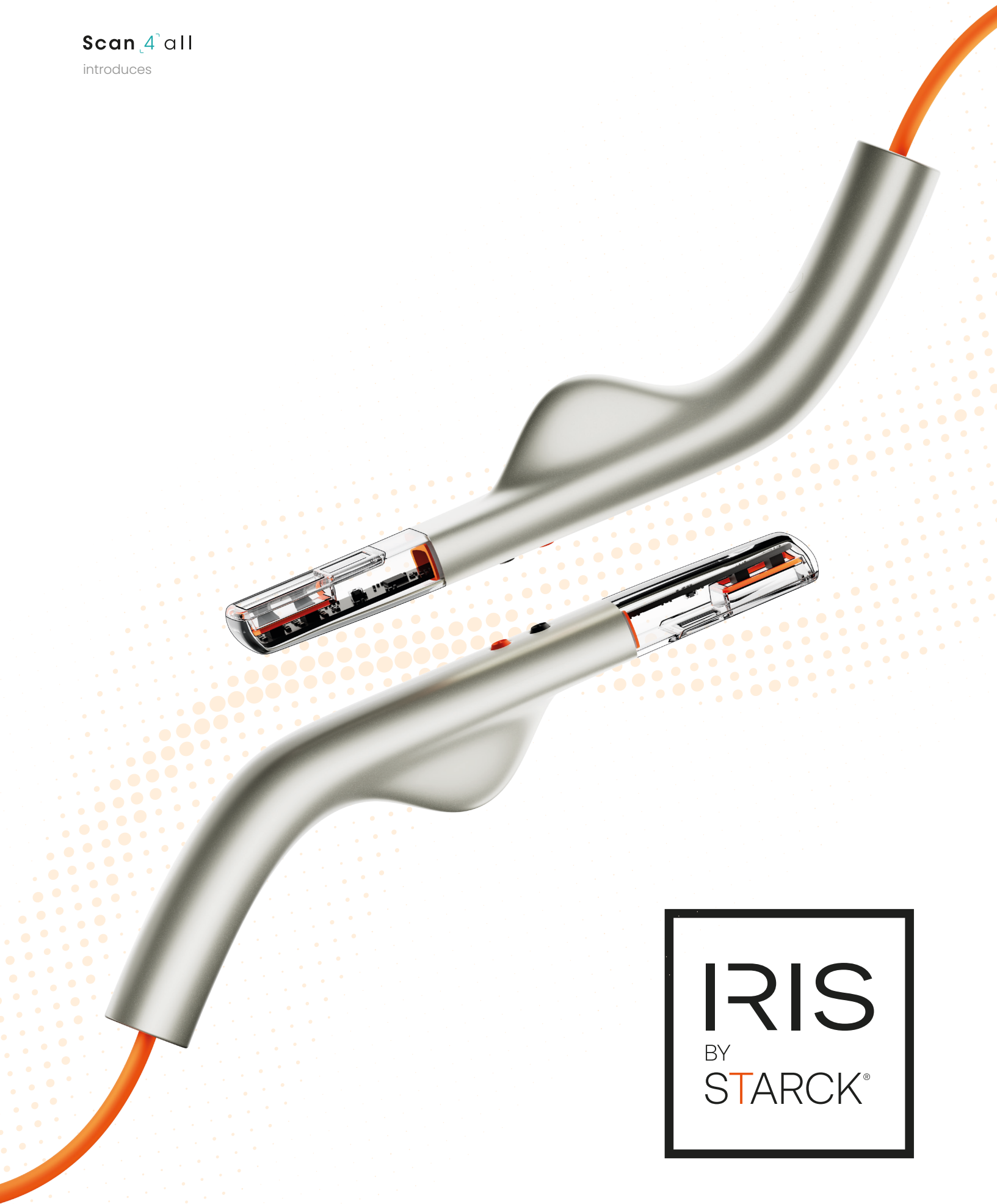

IRIS
BY
STARCK®

INSIDE / **Real time scan**

PREPARATION : THE KEY TO SUCCESS

IRIS BY STARCK® allows you to capture 3D fingerprints in a single step.



The computer must remain plugged into a power outlet during each scan.



Avoid direct light on the area to be scanned as this may interfere with the scanner's lighting calibration and affect its operation.



Dry the mouth and teeth thoroughly for optimal capture quality.




Wear black gloves to minimize unwanted reflections.



Position the handpiece as indicated on the following page. These are your recommended starting points.

USE THE CORRECT MOVEMENT DURING SCANNING

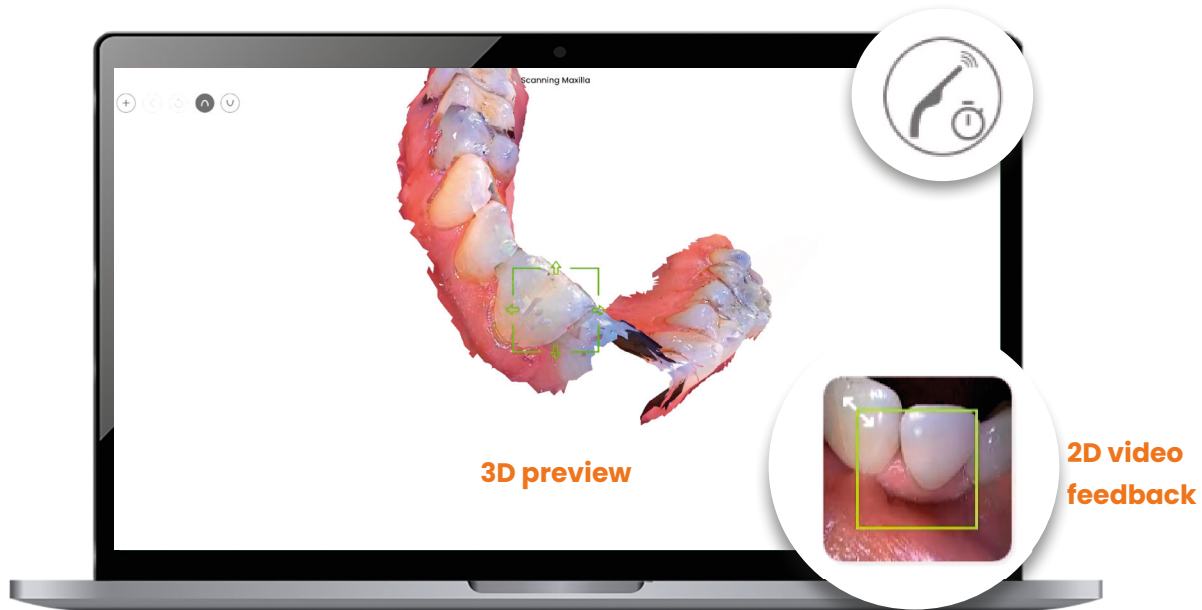

Follow your scan live with the real-time 3D preview.


Scan at a steady pace

Maintain a constant speed and focus on stable movement, **like a paintbrush on a canvas**. Neither too slow nor too fast.



A timer starts at the beginning of the scan and turns yellow, then red if the capture time is too long.


Real-time visual aids

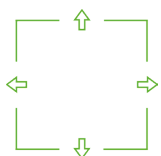

Everything is fine.



Slow down, the scan is stabilizing.

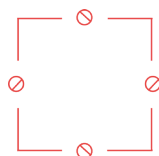

Tracking is lost, return to an area that has already been scanned.

Keep your movements fluid


The IRIS BY STARCK® scanner glides over the teeth smoothly and without rushing. Placing the tip against the teeth ensures the correct distance for an accurate scan.

To begin your scan, you can start with the occlusal surface, then move to the lingual or palatal surface, and finish with the vestibular surface.

For incisors, rock the scanner along the occlusal axis from the inside to the outside and finish by sliding it lightly over the front surface.

SCAN FREELY

Pause at any time to inspect or manipulate your 3D preview and check if the scan is satisfactory. Black areas indicate a lack of data, so remember to complete the scan if necessary.

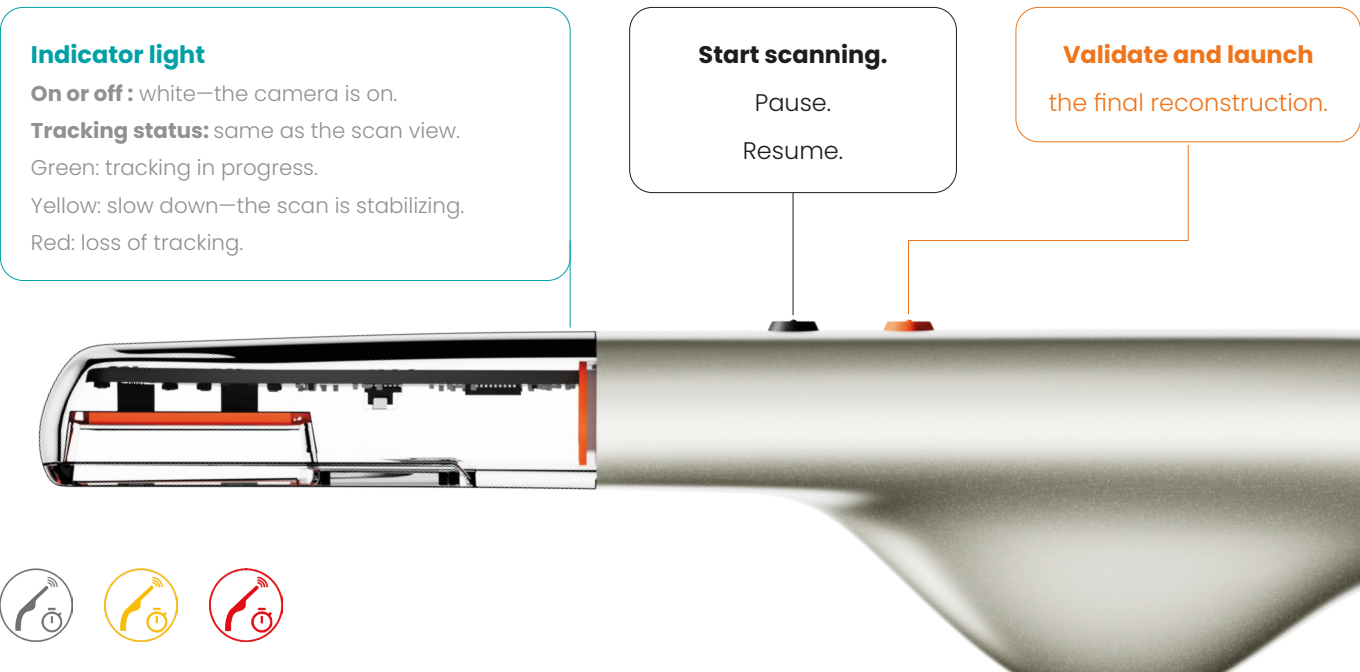

To ensure a smooth scanning experience: use the stopwatch icon (in the top-right corner of your real-time scan view).


Gray: You can continue.

Yellow: Scanning experience is slowed down.

Red: Prioritize starting the final reconstruction. You can resume scanning using the “Resume Scan” tool.

Use the right tools after scanning

Take advantage of scan management features, tools to support you in your clinic, and comparative analysis and measurement options (more information on the file *INSIDE / Tools related to scanning*).

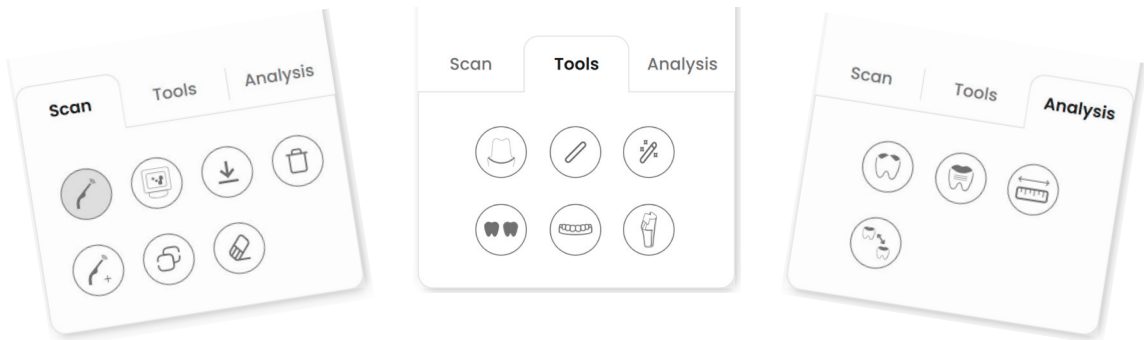

After your scanning process, follow your preferred workflow until the orders are finalized (see the files *INSIDE / My orders ; INSIDE / Open workflow ; INSIDE / Guided workflow*).

GO FURTHER WITH IRIS BY STARCK®.

Explore all the resources available to you.

-  Inside / **Quick start**
-  Inside / **My settings**
-  Inside / **DIRECT SCAN**
-  Inside / **My account**
-  Inside / **System informations**
-  Inside / **Tools related to scanning**
-  Inside / **IRIS Connect**
-  Inside / **SAV**
-  Inside / **My orders**
-  Inside / **Real time scan**
-  Inside / **My notifications**
-  Inside / **Computer specifications**
-  Inside / **My partners**
-  Inside / **Guided workflow**
-  Inside / **My patients**
-  Inside / **Open workflow**



Access all the
documentation.

Manufacturer

Scan⁴all

Excellence, for everyone.

137 rue Garance 31670 Labège - FRANCE

N° SIRET 507 484 863 00059

 www.linkedin.com/company/scan4all/

IRIS BY STARCK®

Class 1 medical device (MD) intended for digital dentistry. CE.

For healthcare professionals only.

© Scan4all - Version 1.0 - 2026