



Scan 4 all
introduces

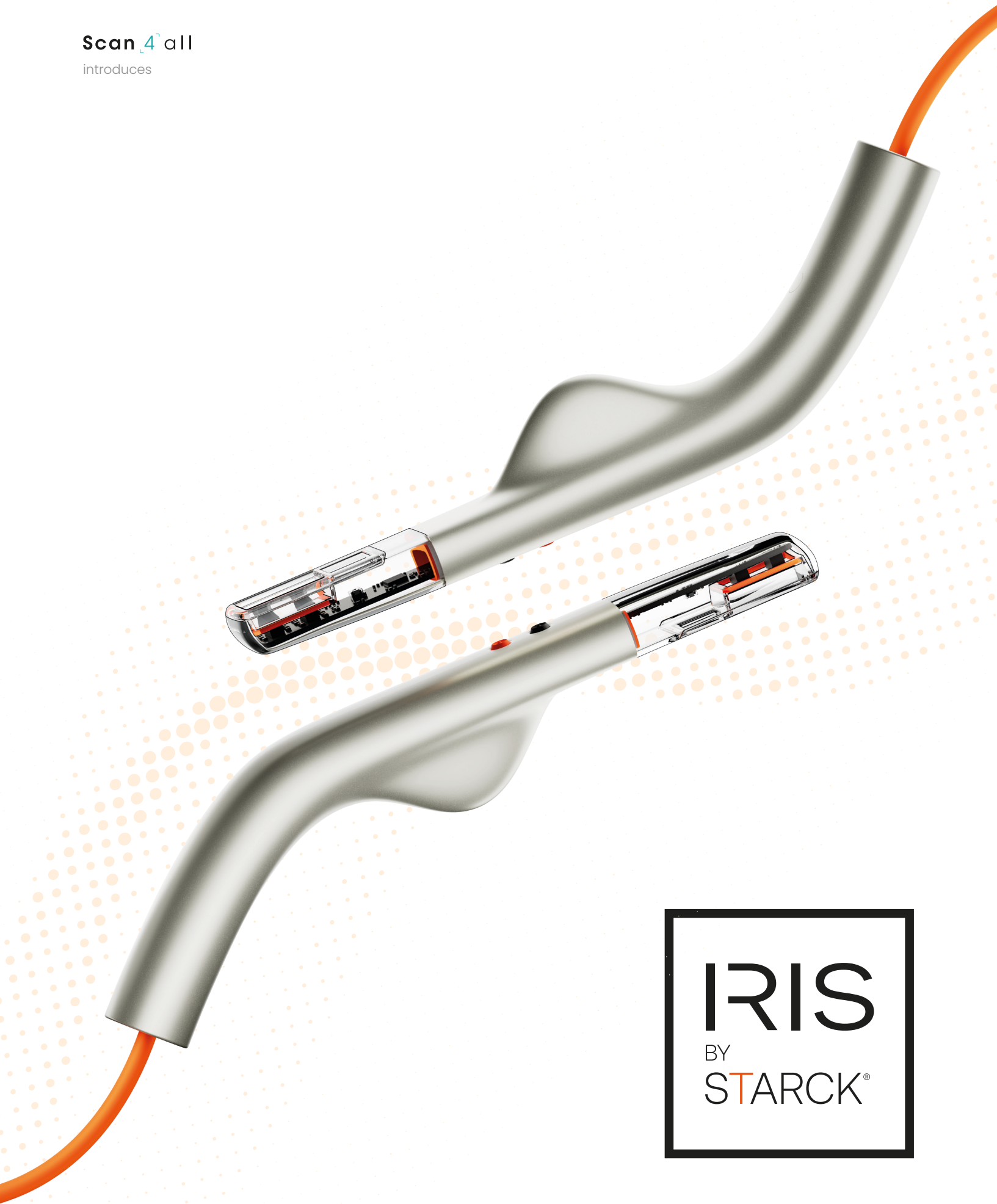

IRIS
BY
STARCK®

INSIDE / Open workflow

THE WORKFLOW, A TAILORED PROCESS

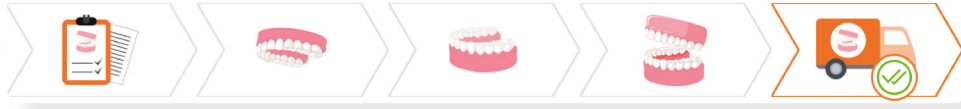


Edit the prosthetic file > **Maxillary** scan > **Mandibular** scan > **Occlusion** scan > Send **order**.

The workflow corresponds to the different steps required to send an order : from creating the prosthetic file to sending the selected scans.

YOUR JOURNEY, YOUR AUTONOMY

In open workflow, each step is independent (when you click on the timeline). You are responsible for performing the scans necessary for a complete order after filling out a prosthetic form.

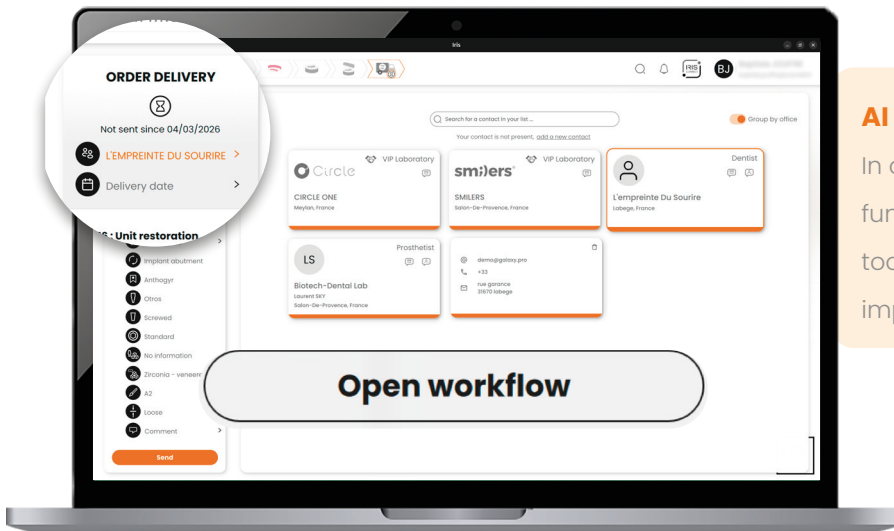

AI implant available

In open workflow, the AI implant function is available via the scan tools and requires you to enter the implants and teeth concerned.

You know which scans you need.

- › You have filled out a prosthetic form and want to select the arches to be scanned yourself.

You already have the necessary impressions.

- › You will place your order directly using the scans from the patient's history.
- › You want to scan an additional arch.
- › You only need an occlusion and want to perform a bite registration.



Why use the open workflow ?

OPTIMIZED FOR COLLABORATION

The DIRECT SCAN function is an open workflow feature that allows you to start a scan immediately, without needing an order in progress.



Maxillary scan > Mandibular scan > Occlusion scan



Maximum time savings


The scans are then attached to the patient file.
The practitioner can complete the prosthetic file and finalize the order later.



Designed for teamwork

Dental assistants can scan independently. Simply select the arch or arches to be scanned.

FINALIZE AN ORDER

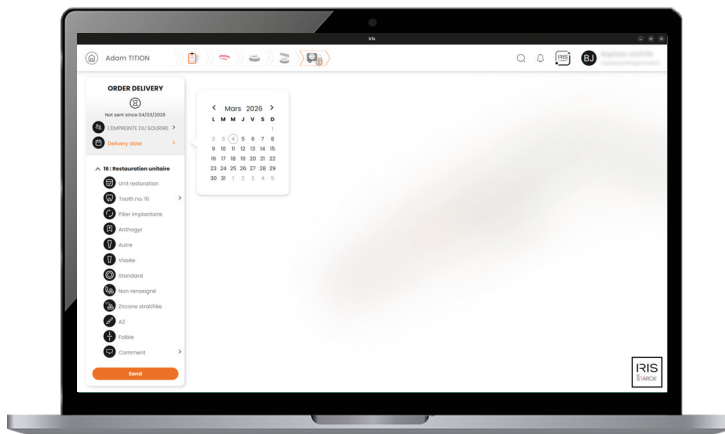

After saving the scan to the patient file, select the patient file and click on the DIRECT SCAN procedure in its history.



Select the laboratory of your choice, complete the prosthetic form and then confirm. Finally, click on the "Order delivery" step and select your delivery date.



IRIS CONNECT keeps you informed.

Dates, exchanges, tracking : everything is centralised and updated instantly. Chat directly with laboratories from the IRIS BY STARCK® CONNECT interface. Smart notifications keep you informed of order progress.

GO FURTHER WITH IRIS BY STARCK®.

Explore all the resources available to you.



Inside / [Quick start](#)



Inside / [My settings](#)



Inside / [DIRECT SCAN](#)



Inside / [My account](#)



Inside / [System informations](#)



Inside / [Tools related to scanning](#)



Inside / [IRIS Connect](#)



Inside / [SAV](#)



Inside / [My orders](#)



Inside / [Real time scan](#)



Inside / [My notifications](#)



Inside / [Computer specifications](#)



Inside / [My partners](#)



Inside / [Guided workflow](#)



Inside / [My patients](#)



Inside / [Open workflow](#)



Access all the
documentation.

Manufacturer

Scan4all

Excellence, for everyone.

137 rue Garance 31670 Labège - FRANCE

N° SIRET 507 484 863 00059

www.linkedin.com/company/scan4all/

IRIS BY STARCK®

Class 1 medical device (MD) intended for digital dentistry. CE.

For healthcare professionals only.

© Scan4all - Version 1.0 - 2026